

# Visit Baby Monitor

BE1491-868



## Features

### Reliable transmission

Up to 250 m operating range

### Flexible

Adjustable sensitivity and delay

### Customizable notification

Use a unique signal for each child

### Optional contact mat

Know if your child leaves the bed

### Easy to use

No complicated menus

**MD**

CE-marked medical device  
under EU MDR 2017/745

## Be there for your baby.

Caring for your baby's health and safety is one of the most important things you do as a parent. Bellman & Symfon offers a fully secure baby monitoring solution for parents with hearing loss, so you are always close to your baby. The Visit Baby Monitor will notify you as soon as your baby needs you, even if you are asleep or a couple of rooms away.

### How it works

When the baby monitor is activated, it signals your Visit receiver that notifies you with sound, flashes or vibrations (depending on the receiver).

### Secure and reliable

Visit uses a secure radio link that covers your entire home, so you always know that you are connected.

### Adjustable sensitivity and delay

You can adjust the sensitivity and delay of the baby monitor so it's not activated when your child moves in her sleep.

### Customizable notification

If you have more than one child you can choose different signals so you know who needs your attention.



# Technical specifications




## Size and weight

- Height: 100 mm, 4.0"  
Width: 65 mm, 2.6"  
Depth: 27 mm, 1.1"  
Weight: 118 g, 4.2 oz (incl. batteries)



## Buttons and connectors

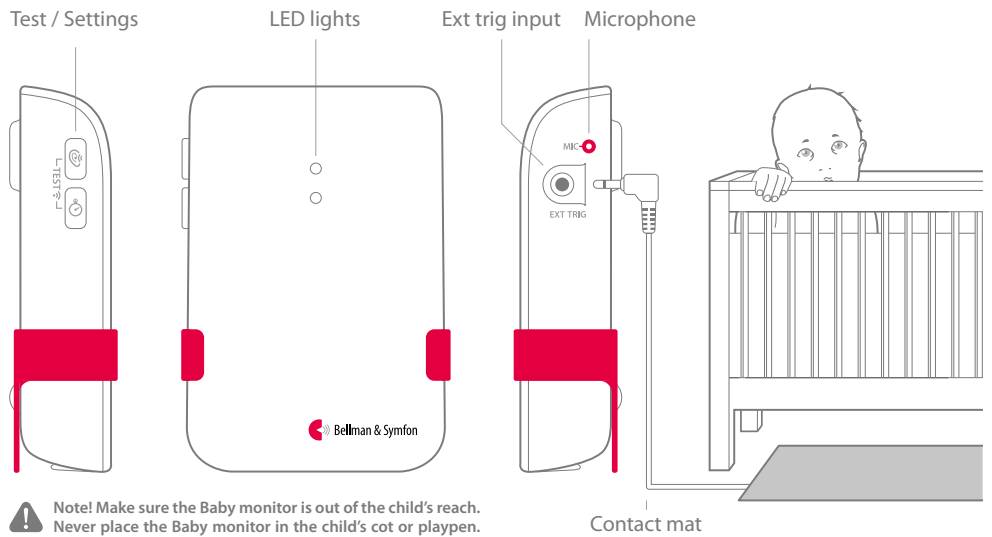

## Inputs/outputs

- 3.5 mm external trigger input

## Activation

- The test buttons and the internal microphone
- The contact mat accessory

## Settings

- Sensitivity: High, Medium, Low
- Delay: 30 sec, 10 sec, 1 sec

## Frequency and coverage

- Radio frequency 868.3 MHz
- Coverage: 50 – 250 m (55 - 273 yd)  
Depending on the characteristics of the building.

## Approvals



This product is a medical device in conformity with Regulation (EU) 2017/745 (Medical Device Regulation).

## Power and battery

- Power supply  
Battery power 2 x 1.5 V AA lithium or alkaline
- Power consumption  
Active < 70 mA, idle position < 400 µA
- Operation time  
Alkaline ~ 6 months, lithium ~ 1 year

## Number of radio keys

- 64 as standard. Up to 256 with optional software.

## Maintenance and cleaning

- Maintenance free. Clean with a dry cloth.
- Do not use household cleaners, aerosol sprays, solvents, alcohol, ammonia, or abrasives.
- Not reusable.

## Environmental requirements

- For indoor use only.
- Temperature: 0 to 35° C, 32 to 95° F
- Relative humidity: 15% - 90% (non condensing)

## In the box

- BE1491-868 Visit Baby Monitor
- Table stand
- 2 x 1.5V AA Alkaline batteries



## Accessories

The following accessories are available:

- BE9026-P02 Contact Mat

© 2025 Bellman & Symfon AB. All rights reserved.

Bellman, Bellman & Symfon and the Bellman & Symfon logo are registered trademarks of Bellman & Symfon AB. This leaflet is subject to change without prior notice.

BE1491-868\_001DSH1.0\_EN

[bellman.com](http://bellman.com)